

DATA SHEET

Product Name: BioPure™ Scramble Amyloid beta 1-42, (**Scramble Aβ42**)
Recombinant, NH₄-Formate salt

Catalog #: AA-0004

Description: This peptide is a scrambled sequence of Amyloid beta 1-42 (AA-0001). Amyloid beta 1-42 (Aβ42) plays a central role in the pathogenesis of Alzheimer's disease (AD). Recombinant Aβ42 has reproducibly been shown to have higher *in vitro* toxicity and aggregate significantly faster than synthetically prepared Aβ42 due to the absence of (n-1) deletion products and racemized amino acids that are characteristic of chemical synthesis^{1,2}.

Sequence: KVKGLIDGDHIGDLVYEFMDSNSAIFREGVGAGHVHVAQVEF

Counter Ion: NH₄-Formate


Source: Recombinant peptide expressed in *E.coli*, processed by proprietary BioPure method; HPLC purified

Peptide Reconstitution: Resuspend lyophilized peptide in 1% NH₄OH or 10 mM NaOH to a concentration of 1 mg/ml. After the material has visibly gone into solution, sonicate for 1 minute to disrupt aggregates. Because the peptide cannot be stored long term in this solution, dilute with 1X PBS or to the desired concentration for assay and/or storage.

Storage: Peptide is shipped at ambient temperature. Store lyophilized or reconstituted peptide at -20°C or lower. Avoid free thaw cycles.

Characterization:

Molecular Weight: 4557.84 Da **% Peak by HPLC:** >99% **Peptide Content:** >90.0%

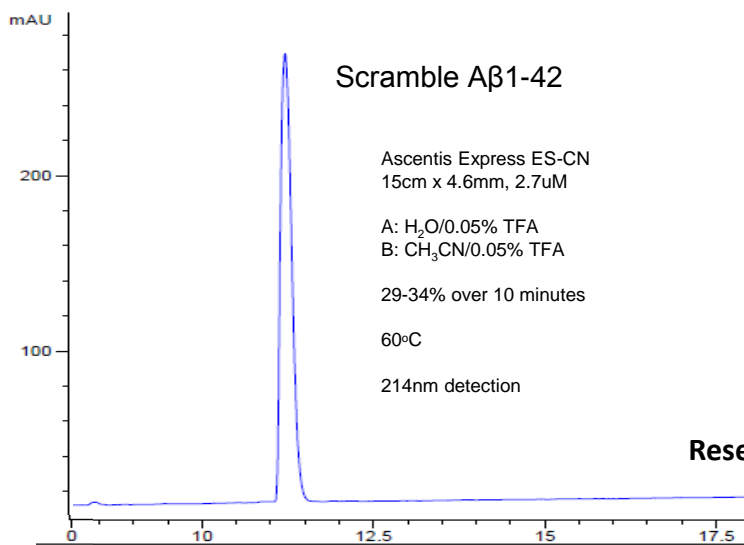

References

1. Bolder SG, Sagis LM, Venema P & van der Linden E (2007) Thioflavin T and birefringence assays to determine the conversion of proteins into fibrils. *Langmuir* 23, 4144–4147.
2. Finder, V.H., Vodopivec, I., Nitsch, R.M., Glockshuber, R. (2007). The Recombinant Amyloid-β Peptide Aβ1–42 Aggregates Faster and Is More Neurotoxic than Synthetic Aβ1–42. *J. Mol. Biol.* 396, 9-18.

Research Use Only – Not for Use in Humans